



Hydrocarbons, C₂ and CO₂

Application Note

Energy & Fuels

Authors

Agilent Technologies, Inc.

Introduction

Gas chromatography with an Agilent PoraPLOT U column and Agilent 490 Micro GC separates carbon dioxide and C₂ hydrocarbons in 25 seconds.



Agilent Technologies


Conditions

Technique : Micro-GC
GC-Channel : Agilent PoraPLOT U, 0.32 mm x 10 m
Temperature : 50 °C
Carrier Gas : Helium, 180 kPa
Injector : 75 ms, 110 °C
BackFlush : after 5 s
Detector : chip TCD

Courtesy : Jim Luong and Rhonda Gras, Dow Chemical Canada

Peak identification

1. carbon dioxide (980 ppm)
2. ethylene
3. ethane
4. acetylene

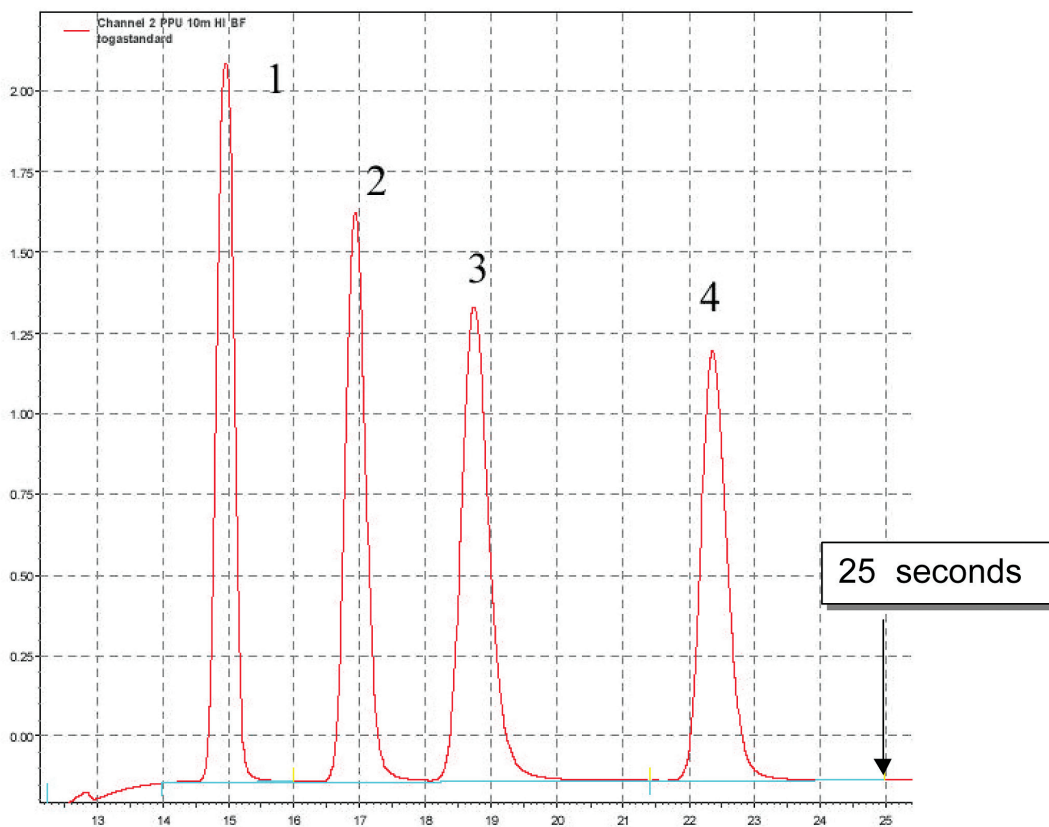

www.agilent.com/chem

This information is subject to change without notice.

© Agilent Technologies, Inc. 2011

Printed in the USA

31 October, 2011

First published prior to 11 May, 2010

A02065



Agilent Technologies